STORZ

KARL STORZ—ENDOSKOPE

Fiberskope

Serien 110xx ..., 111xx ..., 112xx ..., 113xx ...*)

Fiberscopes

Series 110xx ..., 111xx ..., 112xx ..., 113xx ...*)

Fibroscopios

Series 110xx ..., 111xx ..., 112xx ..., 113xx ...*)

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Instruction Manual

Manual de instrucciones

Fiberskope

Serien 110xx ..., 111xx ..., 112xx ..., 113xx ...*)

Fiberscopes

Series 110xx ..., 111xx ..., 112xx ..., 113xx ...*)

Fibroscopios

Series 110xx ..., 111xx ..., 112xx ..., 113xx ...*)

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Important information for users of KARL STORZ instruments

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Thank you for your expression of confidence in the KARL STORZ brand name. Like all of our other products, these products are the result of years of experience and great care in manufacture. You and your organization have decided in favor of modern, high-quality items of equipment from KARL STORZ.

This manual is intended to serve as an aid in the proper use, cleaning, disinfection and sterilization of fiberoscopes manufactured by KARL STORZ. The fiberoscopes are sensitive precision instruments which require particular care during preparation. To prevent damage during preparation, support is provided here for disassembly, cleaning, sterilization and assembly of the KARL STORZ fiberoscopes. For this reason, read these instructions carefully and keep them as a ready reference.

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DCI®-Fiberskope (Direct-Coupled-Interface)

Information regarding this instruction manual

This instruction manual describes three fiberoscope versions, distinguished from one another by the imaging system. If a section refers exclusively to one of these versions, it is thus appropriately labeled.

Telescope

Telescope with special filters for PDD/AF

DCI® fiberoscopes (Direct-Coupled-Interface)

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
KARL STORZ is constantly working on the further development of all products. Please appreciate that changes to the scope for supply in form, equipment and technology are possible for this reason. Therefore, no claims may be derived from the information, figures and descriptions in this manual.

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Controls, connectors, accessories and their uses

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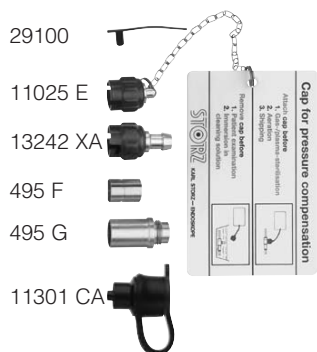- 495 F Light cable adaptor for Wolf
495 G Light cable adaptor for KARL STORZ, Olympus, Winter & Ibe
11025 E Pressure compensation cap (red)
11301 CA Lip valve, disposable version, pack of 20
13242 XA Adaptor for leakage test in cleaning and disinfection machines (green)
29100 Plug for instrument channel, pack of 10, reusable

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
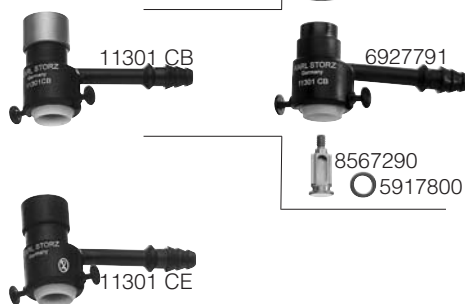
- 11301 D1 LED battery-powered light source with fine screw thread
11301 D2 LED battery-powered light source with click connection
11301 D3 LED battery-powered source with quick-release thread

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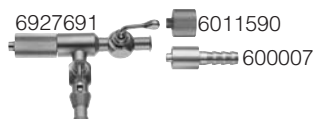- 6927691 Three-way double stopcock adaptor with LUER-lock
6011590 LUER cap
600007 Barbed tube, thin, with LUER-lock

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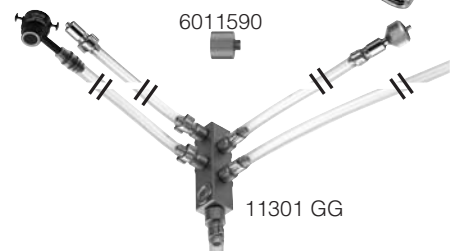
13272 Mouth piece

11301 CF LIPP tube holder

11301 BM Adaptor for leakage test in Belimed cleaning machines

11301 CD Adaptor for machine cleaning the working channel (green)

11301 HH Adaptor for leakage test in BHT cleaning machines

11301 FF2 Adaptor for leakage test in ETD2 cleaning machines with nonreturn valve

11301 GG Adaptor for cleaning the working channel in ETD2 machines

11301 KK Adaptor for cleaning the working channel in ETD3 cleaning machines

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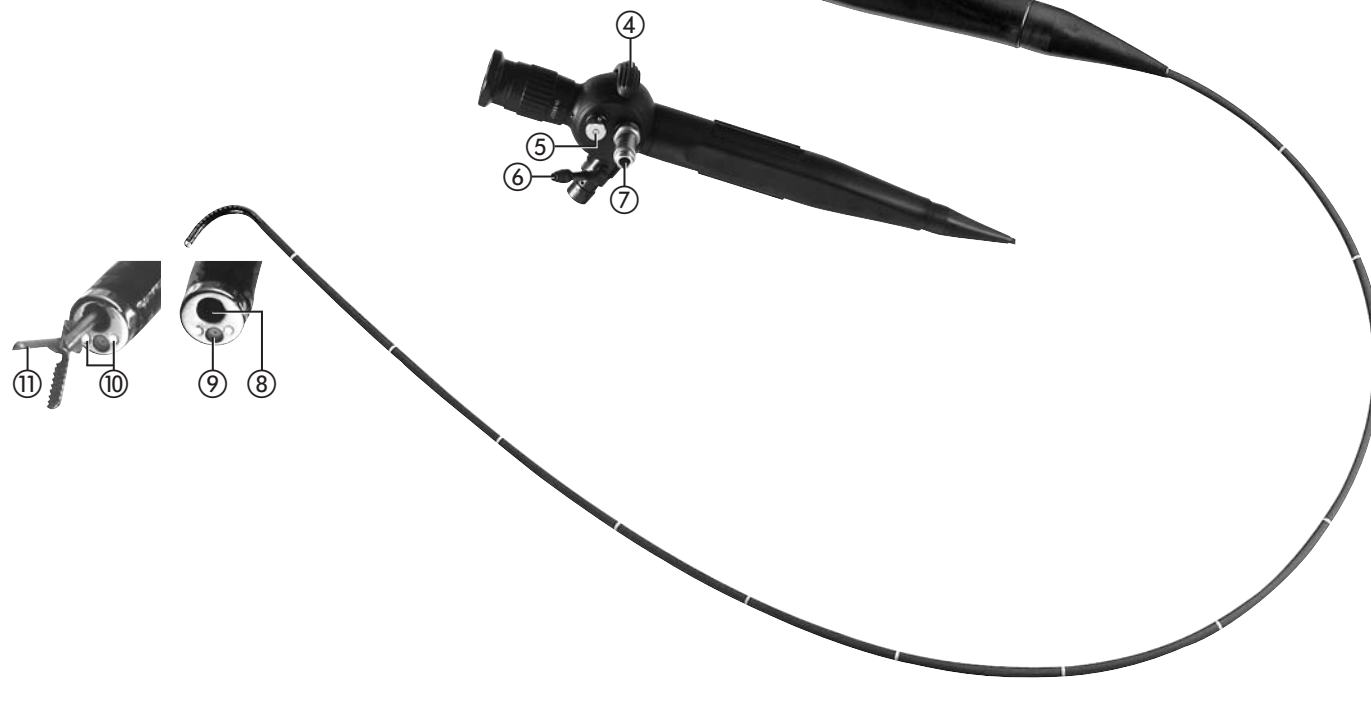
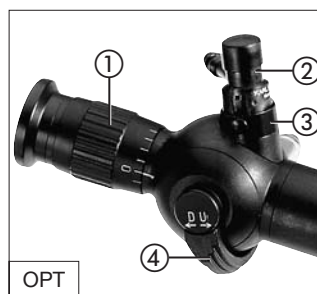
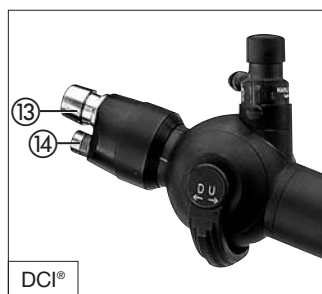
220260030 DCI® adaptor, for connecting endoscopes with DIN eyepiece to DCI® camera heads

20260031 DCI® eyepiece adaptor for connecting DCI® telescopes to camera heads with mounts for DIN eyepieces

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- ① Focusing ring (if available)
- ② Suction valve insert (if available)
Press = suction activated
Release = suction deactivated
- ③ Valve nozzle (if available)
- ④ Up/down control
- ⑤ Valve for pressure compensation
and leakage tester connection
- ⑥ Suction outlet
- ⑦ Light inlet piece
- ⑧ Instrument channel
- ⑨ Objective lens
- ⑩ Light outlet
- ⑪ Instrument in instrument channel

- ⑫ PDD filter changer:
ALA for special tumor marker
based on aminolevulinic acid
W for white light
AF autofluorescence mode

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**Important information for users of
KARL STORZ instruments II**

Fiberscopes Series.....II
Information regarding this instruction manual.....II

**Controls, connectors
and their uses III**

Safety instructions 3

Intended use.....3
User qualification.....3
Safety precautions when
using the fiberscopes.....3
Warnings and cautions.....4

Handling 7

Unpacking7
Basic equipment.....7
Accessories7
Initial operation8
Fitting the suction valve8
Fitting the tube holder.....8
Fitting the lip valve.....8
Handling the fiberscopes9
Using the fiberscopes10

Operating instructions 14

Preparation of the fiberscope.....14
Controlling the distal tip.....15
Suction16
Insufflation.....16
Focusing16
Instrument insertion17
PDD version17
DCI® fiberscope with camera18
DIN eyepiece on DCI® fiberscope.....19
DCI® camera on DIN eyepiece20

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Cleaning, disinfection, care and sterilization..... 21

Preparation Overview	24
Preliminary cleaning	27
Disassembling suction valve insert.....	29
Manual preparation	30
Machine preparation	34
Cleaning connection in EDT3	
cleaning machines.....	36
Cleaning connection in EDT2	
cleaning machines.....	39
Inspection after preparation	42
Assembling suction valve insert.....	43
Gas sterilization/plasma sterilization	44
Chemical sterilization.....	46
Storage	46
Accessories	48
Flexible forceps	48
Papillotome/polypectomy snare	49
Flexible injection cannula	50

Maintenance 51

Servicing and repair	51
Limitation of liability	51
Warranty.....	51
Directive compliance	51
Repair and replacement program.....	52
Important information	52

Appendix..... 53

Chemicals for the preparation.....	53
Lubrication gels	55

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Please read these safety instructions carefully before using the instrument. Make certain you are well acquainted with the function and operation of the instrument before initial use on the patient.

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Intended Use

Fiberscopes are used in diagnostic and surgical endoscopy. They must only be used by appropriately trained and qualified technical staff.

Indications: Examination, diagnosis and/or therapy in conjunction with endoscopic accessories in various medical fields.

Contraindications: No contraindications are known which directly relate to the product. The fibroscope must not be used if, in the opinion of a qualified physician, such an application would endanger the patient, e.g. due to the patient's general condition, or if the endoscopic method as such is contraindicated.

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User qualification

KARL STORZ fiberscopes may only be used by persons with an appropriate medical qualification and who are acquainted with the endoscopic technique.

The information given in these instructions only serves to instruct in the correct handling, cleaning, disinfection and sterilization of the fiberscopes. It is not intended as an introduction to the endoscopic technique.

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Safety precautions when using the fiberscopes

Fiberscopes must be used according to the medical rules and procedures of endoscopy recognized for this method.

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Warnings and cautions

Please read this manual and follow its instructions carefully. The words, **Warning**, **Caution** and **Note** convey special meanings. Wherever they are used in this manual, the text should be carefully reviewed to ensure the safe and effective operation of this product. To make the words **Warning**, **Caution** and **Note** stand out more clearly, they are accompanied by a pictogram.

Warning: A Warning indicates that the personal safety of the patient or physician may be involved. Disregarding a Warning could result in injury to the patient or physician.

Caution: A Caution indicates that particular service procedures or safety precautions must be followed to avoid possible damage to the product.

Note: A Note indicates special information about operating the product, or clarifies important information.

Warning: Before using the fiberscope, read the following safety instructions carefully, to avoid putting your patients, personnel or yourself at risk.

Warning: Inspect the fiberscope prior to each surgical procedure to ensure that it is functioning correctly and that it has been correctly cleaned, disinfected and (where applicable) sterilized.

Warning: This equipment may not be used during the discharge of a cardiac defibrillator. Remove all endoscopic equipment from patient contact and from the surgical site.

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Warning: Verify the sound condition of the fibroscope as well as of accessories used in combination before every application. Damaged fiberscopes or damaged accessories must not be used.

Warning: Combinations of medical devices are only then assured to be safe if

- they are identified as such in the respective instruction manual **or**
- the intended purpose and the interface specifications of the devices used in combination permit this.

Pay careful attention to the instructions and interface specifications of medical products used in combination.

Warning: Lengthy direct contact between the endoscope tip and surrounding tissue may cause tissue damage due to the potentially high temperatures on the tip of the endoscope and must thus be avoided.

Warning: The fibroscope is delivered **unsterilized**, and must therefore be cleaned, disinfected and (where applicable) sterilized prior to initial use and each subsequent reuse.

Warning: Remove the red pressure compensation cap from the vent port

- for every endoscopic intervention
- before immersion in liquid/chemical sterilization.

Warning: Before the fibroscope is removed from the patient, the locking mechanism ⑮ must be released and the distal tip of the fibroscope put in the neutral, non-angled position. Otherwise the instrument may be damaged or the patient injured.

Warning: If a malfunction should occur during application on a patient, stop the application immediately. Place the distal tip of the fibroscope in the neutral, non-angled position and remove the fibroscope slowly and carefully from the patient.

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Caution: Do not steam sterilize (autoclave) the fiberscopes.

Caution: The maximum permissible temperature for cleaning, disinfection, sterilization, transport and storage is 65°C.

Caution: Never immerse fiberscopes in normal saline solution, since even only brief contact can lead to pitting and corrosion.

Caution: Material damage: if the outer casing of the fiberscope is cleaned with alcohol, this will damage the material.

Caution: Fiberscopes must not be cleaned in an ultrasound bath under any circumstances.

Caution: Place the red pressure compensation cap on the vent port before

- gas sterilization/plasma sterilization
- shipping

Caution: Only use KARL STORZ accessories (guide probes, cleaning brushes etc.). Accessories from other manufacturers may damage the fiberscopes.

Caution: Only use guide probes with flexible ends. Other types of guide probe may cause considerable damage to fiberscopes.

Caution: Make certain to read the 'Cleaning, sterilization and care of KARL STORZ instruments' instructions. Methods of cleaning, disinfection and sterilization are explained in detail there.

Note: Any damage to the fiberscope resulting from incorrect operation is not covered by the guarantee.

Note: Follow local governing ordinances and recycling plans regarding disposal or recycling of device components.

Note: If the instrument is in the working channel, its additional rigidity may limit angling.

Note: Separate instruction manuals are available for **flexible videoscopes**; they are therefore not described in this brochure.

Note: Federal law restricts this device to sale by or on the order of a physician (USA only).

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Unpacking

Check for missing items and evidence of shipping damage.
File any complaints to KARL STORZ or the supplier immediately.

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Basic equipment

1 Fiberscope	Art. no.
1 Suction valve insert, disposable version, pack of 20 inserts	11301 CE
1 Leakage tester	13242 XL
1 Pressure compensation cap	11025 E
1 Plug for instrument channel, pack of 10	29100
2 Adaptor for cleaning the working channel in cleaning machines (green)	11301 CD
1 Cleaning brush for 2.8 x 50 and for 2.8 x 70	27651 A*)
for 3.7 x 65	27651 AL*)
for 5.0 x 65	27651 B
for 5.6 x 65	27651 C
for 6.4 x 65	27651 C
1 Biopsy forceps (only bronchofiberscopes)	11001 KL
1 Grasping forceps (only bronchofiberscopes)	11002 KS
1 Three-way double stopcock adaptor with LUER-lock for O ₂ connection (only intubation fiberscopes)	6927691
1 LIPP tube holder (only intubation fiberscopes)	11301 CF
1 Mouth piece (only gastro-, intubation and bronchofiberscopes)	13272

*) For the fiberscopes 11278 A/AI/AU/AIU use the cleaning brush 27651 A (27651 AL is not suitable for these fiberscopes).

Accessories

1 Suction valve insert, reusable version	11301 CB
1 Battery light source	11301 D
1 Adaptor for leakage test in cleaning and disinfection machines (green)	13242 XA
4 Cleaning adaptors for cleaning machines	11301 GG/FF/HH/KK

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Initial operation

Warning: The fiberscope is delivered **unsterilized**, and must therefore be cleaned, disinfected and (where applicable) sterilized **prior to** initial use and each subsequent reuse.

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Fitting the suction valve

Place the suction valve insert ② on the valve nozzle ③ and turn until the projections on the side of the insert lie over the nozzle recesses. Push gently until the suction valve insert clicks into the valve nozzle.

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Fitting the tube holder (only intubation fiberscopes)

Push the tube holder onto the fiberscope from distal and secure it by turning the two parts in opposite directions.

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Fitting the lip valve (only bronchofiberscopes)

Place the disposable lip valve onto the instrument channel ⑧ and tighten gently by turning it in a clockwise direction (plastic part).

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Handling the fiberscopes

The optical section of KARL STORZ fiberscopes consists of a fiberoptic bundle in a flexible and watertight casing, a light guide, as well as objective and eyepiece lenses. Like all optical systems, these parts are very delicate. For this reason, careful handling of the fiberscopes is essential for a long service life and a consistent quality. Therefore be aware of the following points:

- Never kink or tightly coil the fiberscopes and light guides, so that no optical fibers break off. It is therefore best to suspend the fiberscopes for storage or keep in large containers.
- Do not jerk the up/down control ④ or move forcefully.
- Do not insert the instruments or cleaning brushes forcefully through the channels. This may cause damage to the casing, especially if the tip of the instrument is deflected.
- Carefully protect the objective lens ⑨ at the distal end from damage. Do not hit. Never remove soiling with sharp-edged objects.
- Never store or disinfect the fiberscope together with sharp-edged objects to prevent damage to the exterior.
- Never place heavy objects on top of the fiberscope.

Note: Every fiberscope supplied by KARL STORZ includes a photo showing the condition of the fiber bundle upon delivery. A copy of this photo is kept on file in our quality control office. The wedge-shaped marking (at 12 o'clock) serves as a guide for the user.

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Using the fiberscopes

The following instructions are recommended for the proper operation of the fiberscope and accessories. This is not a reference for endoscopic techniques.

Consult the appropriate medical literature for information on endoscopic techniques.

You can find accessories for your fiberscope (catheters, flexible scissors and forceps, guide wires, stone baskets, lithotripter probes, laser fibers, electrodes etc.) in the appropriate KARL STORZ catalog.

Warning: Always adjust the light source to the minimum illumination intensity necessary to achieve optimum visualization of the operating field. The higher the intensity setting of the light source, the greater the heat energy that will be generated at the tip of the fiberscope.

Warning: Never lay the end of the fiber optic cable or the distal tip of the fiberscope on or under a sterile cover cloth. The light intensity may lead to burns on the patient and/or sterile cover cloth.

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When introducing the fiberoscope into the patient, the distal tip of the fiberoscope must be in a straight, nondeflected position (control lock-in can be felt and heard).

The fiberoscope must be handled according to the recognized medical rules and procedures for endoscopic methods.

- If no instrument is inserted, close the instrument channel proximally with a plug (Art. no. 29100).
- If an instrument is to be inserted, position the disposable lip valve (Art. no. 11301 CA) to seal off the instrument channel proximally.
- Instruments especially designed for this may be inserted through the instrument channel of the fiberoscopes. Make certain the accessories have a suitable length and diameter. Unsuitable accessories may damage the fiberoscopes. For questions on accessories, please contact KARL STORZ directly or an authorized dealer.

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Note: When the lip valve is first used, the instrument is subjected to increased resistance. This resistance will decrease the more often the lip valve is used.

- When inserting an instrument into the instrument channel, the fiberoscope must be in a straight, nondeflected position.

Caution: When inserting instruments, increasing resistance is felt with increasing deflection. Forced introduction of instruments can lead to damage to the fiberoscope or the instrument channel.

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Caution: If undue resistance is encountered when inserting an instrument, do not attempt to force the instrument. Withdraw the instrument and apply a water soluble lubricant (See list of approved lubricants in Appendix) to the tip of the instrument and attempt to insert again. Prior to removal of the instrument from the patient, the distal tip must be in a straight position.

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- When using forceps, scissors or other accessories that can be opened or closed, ensure that the instrument is completely closed prior to insertion into the instrument channel. Do not attempt to open the jaws until the distal tip is completely free of the channel. Opening the jaws inside the instrument channel can lead to damage to the fibroscope.

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Warning: When using HF electrodes or laser probes, make sure that the active parts are clear of the instrument channel and in full view of the operating surgeon before power is applied. Activation within the instrument channel may damage the fibroscope and/or injure the patient.

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Warning: Make certain to read the instructions for the laser unit and/or HF generator.

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Caution: If accessories with sharp tips, such as EHL probes or laser probes, are inserted into the instrument channel of a deflected fibroscope, this can damage the fibroscope.

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Caution: When using a stone basket, the stone must not be retracted through the instrument channel but instead the entire fibroscope must be removed. Once the stone is removed, the stone basket may then be taken out of the instrument channel.

Caution: The fiberscopes must be cleaned and then prepared immediately after the intervention in order to prevent incrustations. Dried-on incrustations may jeopardize successful decontamination and disinfection

- A leakage test must be carried out after every intervention, to prevent damage during preparation.

Caution: Use only KARL STORZ accessories (guide probes, cleaning brushes etc.). Accessories from other manufacturers may damage the fiberscopes.

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Note: If the instrument is in the working channel, its additional rigidity may limit maneuverability.

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Preparation of the fiberscope

Before using the fiberscope on a patient, exterior damage and operational problems must be eliminated. Check the outer sheath of the fiberscope. Be especially aware of visible damage to the casing, particularly in the area of the fiberscope tip.

Look through the eyepiece and make certain that you have a clear image. The grid (Moiré effect) should be easily recognizable. If the image is not clear but cloudy, clean the eyepiece and lens with a cotton applicator soaked in alcohol, maintaining hygienic conditions at all times.

DCI® fiberscope image quality can be controlled on screen. If necessary, clean image waveguides, light cables and eyepieces as described above.

If it is still cloudy, it is possible that the objective lens is scratched or that there is moisture inside the telescope. The fiberscope must not be used in this case, but instead sent in for repairs.

Caution: Before every application on the patient, check the instrument for leaks.

Check the fiberscope for leaks before every application on a patient.

Make certain to read the specifications on page 23 of this instruction manual.

Connect these supply lines:

- If necessary, remove the adaptor 495 G or 495 F from the light inlet piece and connect the fiberoptic light cable or battery-powered light source to the light inlet piece ⑦.

Note: Not necessary for battery-powered light sources with coarse threads (11301 D3).

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- suction to suction outlet ⑥
- or O₂ supply (under anesthesia) to instrument channel ⑧ via the appropriate adaptor

Caution: The instruction manuals of the supply units used must be observed.

Switch on the respective supply units.

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Insert the DCI® fiberscope into the DCI® camera head until it locks into place.

Insert the light cable (Art. no. 495 DV) into the DCI® camera head as far as it will go.

Warning: The instruction manuals of the supply units used must be observed.

Connect the light cable and the video signal cable up to the appropriate supply units.

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Controlling the distal tip

The up/down movement of the distal tip on the fiberscope can be adjusted at the control ④: move smoothly and slowly.

Some fiberscopes have:

- contrapositive deflection direction (urology);
- a locking mechanism ⑮ for securing the fiberscope tip (gynecology).

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Models with locking mechanism ⑮: To secure the deflection position of the distal tip, turn the locking mechanism ⑮ slowly to distal until resistance can be felt.

Caution: Forcing the locking mechanism ⑮ further in the same direction may damage the locking mechanism ⑮.

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Warning: Before the fiberoptic instrument is removed from the patient, the locking mechanism ⑮ must be released and the distal tip of the fiberoptic instrument put in the neutral, non-deflected position. Otherwise the instrument may be damaged or the patient injured.

To release the locking mechanism ⑮ from the secure position, turn it slowly to proximal.

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Suction

Press the suction valve ② to start suction (when the suction pump has been connected and switched on).

Caution: There should be no instrument in the suction channel during suction. The instrument channel inlet should be closed so that no secondary air is drawn in.

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Insufflation

Under anesthesia, connect the O₂ hose to the instrument channel:

- using a barbed tube (Art. no. 600007) on the instrument channel or
- if you are using the three-way double stopcock adapter with LUER-lock (Art. no. 6927691), attach the hose to the tube holder on the angled part of the adapter. This ensures that the straight-line passage of the adapter remains free for instrument insertion.

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Focusing

The image is focused with the focusing ring ①.

Note: When connecting photo or video systems recommended by KARL STORZ, the fiberoptic focusing ring must be at the zero position (mark).

DCI® fiberoptic instruments are focused using the adjustment wheel of the DCI® camera (see page 18).

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Instrument insertion

Flexible instruments can be inserted through the instrument channel (8) as required.

If the instrument channel is not being used: replace the sealing cap.

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PDD

PDD version

Fiberscopes intended for use in Photodynamic Diagnosis (PDD) or autofluorescence (AF) are illuminated with a special diagnosis light (D-Light C). Adjust the PDD filter changer (12) for special filters designed for this light and possibly for marker substances:

- ALA for a special tumor marker based on aminolevulinic acid
- W for white light
- AF autofluorescence mode

Note: If the filter changer is set to 'ALA' or 'AF', a yellowish image is produced when viewing through the integrated special filter.

Note: No filter changer is available on PDD fiberscope art. no. 11272 CI. The filter (ALA) is permanently installed.

For optimal illumination, use the fluid light cable 495 FP (2.5 m) or 495 FO (1.8 m) between the light source and PDD fiberscope.

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DCI® fiberscope with camera

Assembly

Slide the DCI® into the snap-in DCI® coupling mechanism of the DCI® camera until the coupling snaps in automatically.

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Connecting the special light cable 495 DC or 459 DV to the DCI® camera head: Insert the light cable plug into the DCI® camera head until it locks into place.

While doing this, hold the light cable by the handle piece. Never pull on the cable!

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
Note: The entrance and exit faces of the light cable and endoscope must be kept clean if optimal light delivery is to be maintained. Clean the surfaces with a cotton applicator soaked in cleaning solution, then disinfect it with a cotton applicator soaked in 70% isopropyl alcohol and dry.

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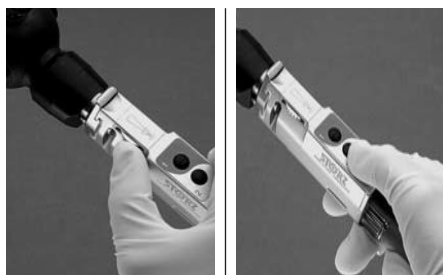DCI® fiberscopes are focused using the adjustment wheel on the DCI® camera.


The camera control unit is controlled using the two programmable control keys (silver, blue) on the DCI® camera head.

Note: Further controls are described in the instruction manual of the camera control unit used.

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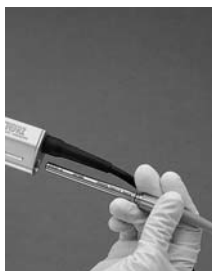Disassembly

Hold the DCI® camera and the DCI® fiberscope still. Push the release button on the DCI® camera: the „snap-in“ DCI® coupling mechanism is released.

Pull the light cable out of the DCI® camera head.

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DIN eyepiece on DCI® fiberscope

The DCI® eyepiece adaptor 20260031 ensures the connection of DCI® telescopes to camera heads with DIN eyepiece mounts.

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Assembly

Slide the DCI® fiberscope into the „snap-in“ DCI® coupling mechanism of the DCI® eyepiece adaptor until the coupling snaps in automatically.

Connect a light cable or a battery-powered light source.

Warning: The instruction manuals of the supply units used must be observed.

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DCI®

Disassembly

Remove the light cable or battery-powered light source.

Hold the DCI® eyepiece adaptor and the DCI® fiberscope. Press the release button on the DCI® eyepiece adaptor. The „snap-in“ DCI® coupling mechanism is released.

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DCI® camera on DIN eyepiece

Assembly

Put the focusing ring (if available) of the fibroscope in the zero position (mark).

Slide the DCI® adaptor 20260030 into the 'snap-in' DCI® coupling mechanism of the DCI® camera until the coupling snaps in automatically.

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Connect the standard light cable at the light inlet piece ⑦ of the fibroscope.

Note: The special light cable 495 DC or 495 DV cannot be used and may have to be removed from the DCI® camera.

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Caution: No further tightening of the outer ring to position the eyepiece is required. This may, among other things, cause damage to the instrument coupler!

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OPT

Disassembly

Hold the DCI® camera and the fibroscope still. Open the instrument coupler for the DCI® adaptor and remove the fibroscope.

Hold the DCI® camera and the DCI® adaptor still. Push the release button on the DCI® camera: the 'snap-in' DCI® coupling mechanism is released. Remove the DCI® adaptor.

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Caution: The instructions 'Cleaning, sterilization and Care of KARL STORZ Instruments' must be followed. There, the cleaning, disinfection and sterilization procedures are explained in detail.

Caution: For cleaning and disinfection, use a product from the KARL STORZ approved chemicals list (see Appendix, pp. 53, 54). The current list of approved products can be found on our website (www.karlstorz.com).

Caution: When preparing and using the solutions, follow the manufacturer's instructions paying close attention to proper concentration and exposure times. Prolonged immersion and incorrect concentrations may result in damage. Bear in mind the microbiological range of action of the chemicals used.

Warning: The stated cleaning, disinfection and sterilization methods alone do not guarantee proper preparation or sterility. This can only be achieved if the staff follows recognized and validated preparation methods.

Warning: For all cleaning of contaminated endoscopes and instruments, liquid-proof protective gloves must be worn. When using a cleaning gun, safety goggles are necessary. Aerosols should be removed using a suction device.

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Cleaning, care, disinfection and sterilization

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Warning: Remember that sterilization can only be performed successfully if

- all organic and chemical residue on the instruments has been eliminated. Sterilization is only effective on clean, disinfected surfaces!
- work proceeds according to recognized and validated methods.
- any deviation from recommended procedures must be validated by the user.
- the recommended sterilization parameters are maintained. These are only valid with sterilization equipment that is properly maintained and calibrated.

Caution: Never clean fiberscopes with a metal brush.

Caution: Constantly changing between different preparation methods causes particular wear on the materials and must therefore be avoided. A once-only change in methods however, e.g. when acquiring new cleaning equipment, is perfectly safe.

Caution: Fiberscopes must not be cleaned in an ultrasonic bath.

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Caution: Never immerse fiberscopes in normal saline solution, since even brief contact can lead to pitting and corrosion.

Caution: Final rinsing during manual preparation of fiberscopes should be carried out using microbiologically pure/sterile water.

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Preliminary cleaning

Immediately after use, remove any heavy soiling from the exterior of the fiberoptic (secretions etc.). Use a disposable cloth moistened with disinfectant (preferably alcohol-free). Submerge the distal end in a container with cleaning or disinfectant solution and activate the suction to clean the working channel.

Caution: Before soaking in disinfectant, check the instrument for leaks. Otherwise fluid penetrating the fiberoptic could cause damage.

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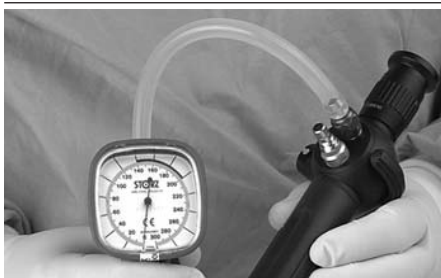
Connect the leakage tester (art. no. 13242 XL) to the valve for pressure compensation ⑤ and build up pressure until the pointer is in the blue marked range. Watch the pointer for approx. 30 seconds: the pressure must not drop during this time. Press the red button on the pressure valve and release the pressure. Remove the test unit from the fiberoptic.

Note: During inflation, make certain that the red pressure relief button is completely closed.

If there is a loss in pressure, the fiberoptic has a leak and must neither be used nor immersed in liquid. In this event, send the fiberoptic in for repairs.

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Warning: Risk of infection: In order to prevent infection, it is strictly forbidden to ship contaminated medical devices! All medical devices must be **decontaminated on-site**, to avoid contact and aerogenous infections among staff. To eliminate the risk of infection, before the fiberscope is dispatched it should be **wiped-down disinfected and sterilized in the NT procedure with the pressure compensation cap in place**. We reserve the right to return contaminated instruments/equipment to the sender.

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- Remove all instruments from the instrument channel ⑧.
- Remove adapter and other removable parts from the light inlet piece ⑦ and instrument channel ⑧.
- Remove suction valve insert ② from the valve nozzle ③: press the projections on the side together and pull out the suction valve insert.
- Dispose of disposable parts removed. If necessary, disassemble other parts further and prepare separately.

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Disassembling suction valve insert

Note: Only the multi-use suction valve insert (silver-coloured cap, metal plunger) is suitable for preparation.

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Manual preparation

Note: Before cleaning, remove both adaptors from the light inlet piece ② to permit free access to the fiber surface.

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Manual cleaning

Submerge the instrument in a cleaning solution approved by KARL STORZ and clean it using a cleaning agent (brush or cleaning gun, Art. no. 27660). Avoid splashing with the cleaning agent.

Caution: Pass the brush through the channel in only one direction and pull back only after the brush head has completely left the distal end. The brush must not be pulled back and forth within the channel, since this can lead to damage.

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Clean the brush head appearing at the distal end with another brush or compressed air before retracting it. Avoid splashing. Then thoroughly rinse the channels repeatedly using cleaning and disinfection agents. If required, add a cleaning agent booster.

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Detach the adaptor, clean the stopcocks and grease (see the 'Cleaning, sterilization and care of KARL STORZ instruments' instructions).

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The adaptor 11301 CD can be used for manual rinsing of the suction/irrigation channel.

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Rinse to remove the cleaning solution.

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Clean the eyepiece and objective lens with a soft cloth, sponge or cotton swab, soaked with 70% alcohol.

Caution: Cleaning with the special cleaning paste should only be carried out if the image as viewed through the endoscope is visibly cloudy or if there are visible deposits on the optical surfaces, it should certainly not be used routinely during preparation.

Stubborn soiling and deposits may be removed with cleaning paste (art. no. 27661). To do this, a small amount of cleaning paste is placed on a clean, moistened cleaning swab and distributed over the glass surface. Through polishing, even stubborn deposits can be removed.

Note: The instrument must be prepared again following treatment with the cleaning paste!

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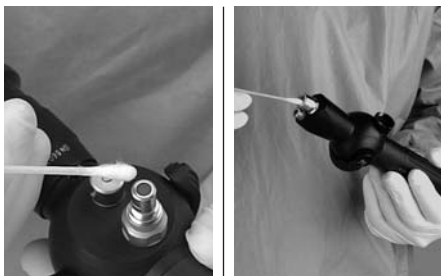
Finally, dry the surfaces of the fibers and the optical end faces with a cotton-tipped applicator soaked in 70% alcohol.

Note: The fiber surface in the light inlet piece ⑦ must be dried off especially thoroughly with alcohol. Residue from disinfectants and detergents in the light inlet piece could bake on if the fiber optic light cable is connected and severely impair light transmission.

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Manual disinfection

All parts of the KARL STORZ fiberoptic scope can be disinfected. For this purpose, materials approved by KARL STORZ and recommended for use with endoscopes should be used. Please refer to the Appendix (Internet).

For disinfection, we recommend a plastic basin (Art. no. 27647 V), in order to avoid scratches on the optical end faces in particular. A plastic basin also prevents electrolytic corrosion, which can occur if different metals are present in a solution.

For disinfection, submerge the fiberoptic scope completely. Ensure that the lumina of the channels and sheaths are completely filled with the solution. There must be no more air (air bubbles) in the channels. Make sure that all the stopcocks, the valve insert and the valve cap are open or have been disassembled.

Immerse the valve insert, adaptors etc. separately.

Rinse out the instrument channel ⑧ with a syringe filled with disinfectant solution. Follow the manufacturer's instructions for concentration and exposure time. Rinse all parts of the instrument through thoroughly at least three times with the prepared solution. The rinsing time should be at least a minute in each case. The rinsing water used must be disposed of since it is contaminated with cleaning solution.

For each rinse, use new microbiologically pure/sterile water. It is absolutely necessary to rinse the fiberoptic scope with microbiologically pure/sterile water in order to prevent potential toxic impact from any remaining chemical residue and to remove surfactants, which could cause the instrument to become conductive. This is also required to optimize the results of sterilization and avoid damage to the surface.

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Dry the outside of the instruments with a lint-free sterile cloth and/or sterile compressed air.

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Remove any water in the instrument and suction channels using sterile compressed air.

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Caution: The pressure must not exceed 50 kPa (0.5 bar/7.25 psi). The cleaning gun might therefore have to be operated with a pressure reducer.

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
Caution: Check the fibroscope for leaks before every preparation in semi- or fully automatic equipment.

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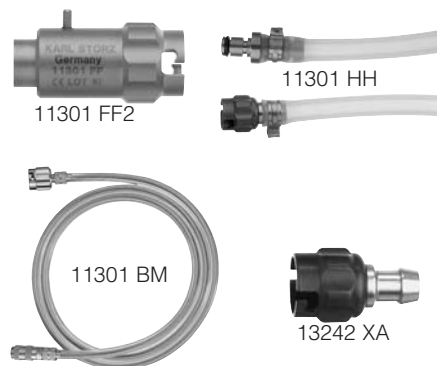Leakage tester connection on washer/disinfector:

- ETD3/ETD2 washers: Adaptor for leakage test (Art. no. 11301 FF2, see pp. 34/37)
- BHT washer/disinfector: Adaptor for leakage test (Art. no. 11301 HH)
- Belimed washer/disinfector: Adaptor for leakage test (Art. no. 11301 BM)
- On the fibroscope using the green adaptor included in delivery (Art. no 13 242 XA); On the machine: consult the manufacturer of the washer/disinfector.

If the preparation machine cannot be connected to the leakage testing unit for any reason, check for leakage manually as described on page 27 or check the adaptor connection using the 13242 XA.

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To prepare the suction/irrigation channel in a washer, connect the green washer adaptor (Art. no. 11301 CD) to a preparation tube and place it on the valve connections.

For preparation in ETD3 machines, use the special adaptor (Art. no. 11301 KK, see p. IV).

For cleaning in ETD2 machines, use the special adaptor (Art. no. 11301 GG, see p. IV).

For fiberscopes with a working channel without a suction valve, connect the washer/disinfector to the working channel ⑧.

To avoid pressure loss in the washer/disinfector, connections which are not in use should be sealed with a LUER-lock (Art. no. 6011590).

Warning: Please observe the instruction manual of the washer/disinfector.

Lay the fiberscope carefully in the tray of the machine. Connect the rinsing tubes to the channels. Proceed with the endoscope cleaning and disinfection program according to the regulations of the machine manufacturer.

Caution: After the program has finished, the fiberscopes are highly susceptible to corrosion due to the ambient conditions in the preparation machine.

Remove the fiberscopes from the machine as soon as possible after the program has finished.

Then dry the optical surfaces with a cotton-tipped applicator soaked in 70% alcohol.

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**Connection for machine preparation
in ETD3 machines**

Dismantling the adaptor 11301 GG:
Press the button corresponding to each tube and
pull the tube off simultaneously.

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1. Slide the free tube end ① of the adaptor (Art.
no. 11301 GG) onto the olive connector of the
green adaptor (Art. no. 11301 CD).

Note: For fiberscopes with a working channel
without a suction valve, pull the tube not
required, e.g. ① out of the metal block.

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4. For machine preparation, attach the green
adaptor (Art. no. 11301 CD) to the valve
connections of the fiberscope.

5. Attach the tube with LUER-lock to the
fiberscope instrument channel ⑧ and secure it.

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6. Insert the other tube ends in the metal part of
the adaptor (Art. no. 11301 GG).

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7. Connect the adaptor (Art. no. 11301 FF) to the 'valve for equalizing pressure and connecting leakage tester'.

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8. Connect the leakage tester to the EDT3 using the adaptor 11301 FF.

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9. Connect the adaptor for cleaning the instrument channel in EDT3 washers (Art. no. 11301 KK) to the connector of the supply tube.

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10. Lay the instrument carefully in the tray of the EDT3 machine.

11. Connect the rinsing tubes to the EDT3 machine.

Warning: To avoid pressure loss in the EDT3 machine, connections which are not in use should be sealed with LUER-locks.

12. Screw the adaptor 11301 KK onto the designated connections of the EDT3 machine, according to the figure.

13. Use the cleaning and disinfection program for fiberoptic according to the EDT3 manufacturer's instructions.

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After cleaning has been completed, remove the green adaptor (11301 CD) for cleaning the working channel by pressing together the projections at the side.

Note: The green adaptor for cleaning the working channel may be left on the tube for later use.

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Remove the leakage tester connection from the adaptor 11301 FF.

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
Caution: The fiberscope could be damaged by liquid which may enter the device via any suction valve inserts (11301 FF) left on the fiberscope after machine preparation. Therefore, you should always dismantle the adaptor (11301 GG) immediately after cleaning is completed!

Pull out the suction valve insert (11302 FF).

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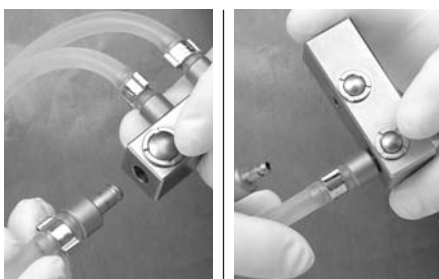To disassemble adaptor 11301 GG: Press respective button and pull off tube simultaneously.

Note: If, according to the high risk classification and/or the Infection Protection Act, sterilization of flexible fiberscopes is not required, the instruments should be stored safely in a hanging position (distal tip pointing down), so that any remaining moisture may drain out.

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**Connection for machine preparation
in ETD2 machines**

Push free end of the adaptor tube (Art. no. 11301 GG) onto the barbed tube (Art. no. 11301 CD).

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Place the adaptor 11301 CD on the valve nozzle ③ of the fiberscope.

Insert the other end of the tube into the metal part of adaptor 11301 GG.

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Place the LUER-lock on the working channel ⑧ of the fiberscope.

Insert the other end of the tube into the metal part of adaptor 11301 GG.

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Insert the supply tube with the fitting for ETD2 cleaning machines into the metal part of adaptor 11301 GG.

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Attach adaptor 11301 FF to the valve ⑤ for pressure compensation and the leakage tester connection of the fiberscope.

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Connect leakage tester connection of the ETD2 cleaning machine.

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Lay the endoscope carefully in the machine basket. Connect the irrigation tubing to the channels. Proceed with the endoscope cleaning program according to the regulations of the machine manufacturer.

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Caution: After the program run, the fiberscopes are greatly at risk of corrosion in the preparation machine.

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Once the program is complete, preferably remove the fiberscopes from the machine immediately.

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Then if necessary dry off the optical surfaces with a cotton swab soaked with 70% alcohol.

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Inspection after preparation

Before assembly it should be checked that all parts and surfaces are undamaged, free from any organic or chemical residue as well as completely clean and dry.

Replace defective parts or send in for repair. Before sending in for repair, all parts must be cleaned and sterilized/disinfectad.

If residue is still present after preparation, prepare the parts again, subject them to special treatment, discard or repair. Remove stubborn deposits on glass surfaces with a special cleaning paste (Art. no. 27661).

The part must then be submitted to a complete preparation cycle.

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
Check plastic parts for discoloration, porosity and flexibility. If the parts have changed substantially from their new condition, they must be replaced. This applies in particular to seals and rubber caps.

Test the instruments for proper operation: check the telescope for a clear view, the instrument for steerability and all canals for patency.

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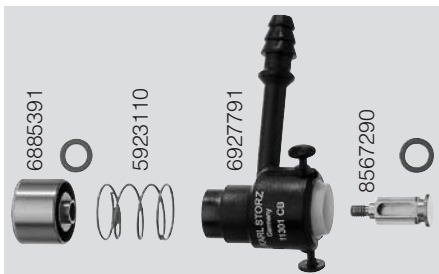Assembling suction valve insert

Note: Only the reusable suction valve insert (silver-coloured cap, metal plunger) is suitable for preparation.

Check seal rings.
Replace damaged seal rings.

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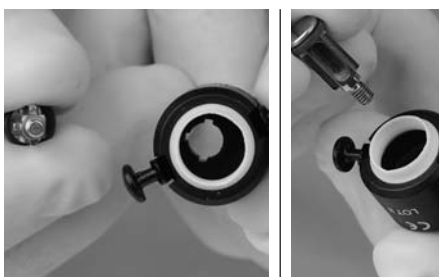
Place plunger in the housing correctly
(use three guide rails).

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Place spring on the upper part of the housing.

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Gas sterilization/plasma sterilization

Sterilization using formaldehyde or ethylene oxide and plasma sterilization are gentle and therefore recommended methods of sterilization.

Note: Observe national laws and regulations for the gas sterilization method.

Note: Fiberscopes which are to be used for interventions in areas of the body which lack microbial colonies, must be sterilized in sterile packaging and stored in a closed cupboard free from contamination and dust.

Caution: If the instrument is to be sterilized, the entire system must be opened by fitting the red cap on the pressure compensation part to ensure pressure compensation.

Warning: Before commencing gas sterilization, the instrument must be completely dry or else sterile conditions will not be achieved.

- Place the instrument securely on a sterilization sieve.

Ethylene oxide gas sterilization has been validated under the following parameters:

Gas mixture	EtO: HCFC*
Temperature (°F, °C)	130°F ± 5°F; 54°C ± 2°C
Relative air humidity	60 ± 20%
Pressure (p _b , kpa)	58 ... 70 (0.56 ... 0.7 bar/8 ... 10 psi)
Exposure time	120 mins
EtO concentration	600 ±30 mg/L

* Gas mixture consisting of 10% ethylene oxide and 90% chlorotetrafluoroethane (HCFC-124; percent by weight)

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For gas sterilization with ethylene oxide, the following airing times must be observed for instruments due to gas absorption of the materials:

	at 20°C	at 42°C	at 55°C
Telescopes	24 hrs	18 hrs	12 hrs
Rubber items	7 days	4 days	3 days
Fiberscopes/ Plastics	5 days	2 days	30 hrs

(minimum time in simple transparent packaging and an application time on the patient of more than 30 minutes).

Maximum acceptable levels (per 21 CFR 812.100; Federal Register, Vol. 43, No. 122, June 23 1978) of residues following EtO sterilization are as follows:

Ethylene oxide:	25 ppm* / 250 ppm**
Ethylene chlorohydrine:	25 ppm* / 250 ppm**
Ethylene glycol:	250 ppm* / 5000 ppm**

** for instruments in contact with skin and mucous membrane,
* for instruments in contact with blood.

The instruments can be used immediately after sterilization with formaldehyde or plasma sterilization. Extra airing time is not necessary.

Remember if you use plasma sterilization that there are limits regarding the size of the lumen and material in regard to instruments which may be sterilized using STERRAD®.

Statements as to the efficiency of sterilization for the respective type of instrument (depending on the length and diameter of the lumina) may only be made by the manufacturer.

Warning: Sterility is only assured if the manufacturer's instructions for the STERRAD® procedure are followed.

Warning: All deviations from the sterilization parameters recommended for the STERRAD® system must be validated by the user.

Caution: Attention must be paid to the information on preparation, sterilization and storage in the "Cleaning, sterilization and care of KARL STORZ instruments" instructions. The methods for cleaning, disinfection and sterilization are explained there in detail.

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Chemical sterilization

Chemical sterilization (STERIS®) is a special method which is particularly suited for thermolabile goods such as fiberscopes. Except for telescopes with a Hamou eyepiece drive mechanism, all KARL STORZ telescopes are suitable for STERIS® sterilization. Color changes to the surface of aluminum alloys could occur during sterilization, which do not however affect the function of the instrument. It must be pointed out in this context that anodized aluminum parts are not always recognizable as such. Sterilization using STERRAD NX is also possible.

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Warning: Sterility is only guaranteed if the manufacturer's instructions for the STERIS System have been followed. This also applies to use of the sterilization containers.

Warning: Chemical sterilization is not possible on greased or oiled surfaces.

Warning: Parts with a lumen must be prepared together with the STERIS Quick Connect kit. All parts must be thoroughly cleaned before sterilization using the STERIS PROCESS™, as also applies to the other sterilization methods.

Caution: Attention must be paid to the information on cleaning, sterilization and storage in the "Cleaning, sterilization and care of KARL STORZ instruments" instructions. The methods for cleaning, disinfection and sterilization are explained there in detail.

Caution: For miniature endoscopes with a removable eyepiece, the eyepiece must either be sealed with a plug or firmly attached to the endoscope before sterilization.

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Storage

Caution: Do not expose fiberscopes to direct sunlight or excessive heat.

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Sterilized fiberscopes should be packaged and stored in sterile cloths or suitable containers in accordance with EN 868, and must be absolutely dry before being put into storage.

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Storage must take place in facilities with normal temperature and air humidity levels.

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
You can find special holders for this in the “Cleaning, sterilization and care of KARL STORZ instruments” instructions.

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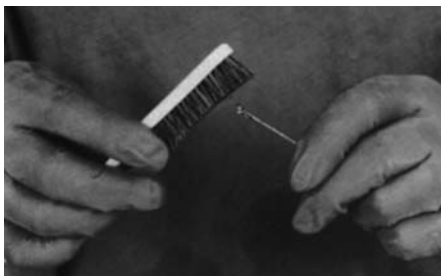Accessories

Flexible forceps

Carefully wipe the outer surfaces of the flexible forceps with a disposable cloth soaked in disinfectant. Clean the forceps jaws with a brush.

Next, clean in an ultrasound bath.

Afterwards, rinse with microbiologically pure/sterile water.

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For handles which can be dismantled:

Dismantle handle as follows:

- Undo the thumbscrew on the handle and remove the rear part of the handle.
- Undo the thumbscrew on the connecting piece and remove the front part of the handle to the rear.

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Also oil the joints of the jaws by holding the instrument so that the oil can flow into the instrument channel and wet the pulling wire. Then re-assemble in reverse order.

The sterilization of flexible forceps is compulsory!

Note: Only use special oils whose lubricating effect remains even in autoclaves (art. no. 27656).

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Papillotome/Polypectomy snare

Dismantle the device completely for preparation.

- Unscrew thumb handle.
- Undo contact socket.
- Unscrew grip piece.
- Remove inner guide sleeve.
- Pull papillotomy or polypectomy loop out of the plastic probe at the front.

Then reassemble the polypectomy snare in reverse order.

Wipe the exterior of the instruments with a cloth soaked in disinfectant. Clean the forceps jaws with a brush. Then clean in an ultrasound bath. Afterwards, rinse with microbiologically pure/sterile water.

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
Assembly of papillotomes:

- Screw inner guide sleeve into the grip piece
- Fit contact socket.
- Insert loop into plastic probe as follows: insert long end of loop into the rear opening of the probe and push to the rear, keeping the probe stretched.
- Insert shorter end of loop into the front opening of the probe and slide both parts towards the rear keeping them parallel.

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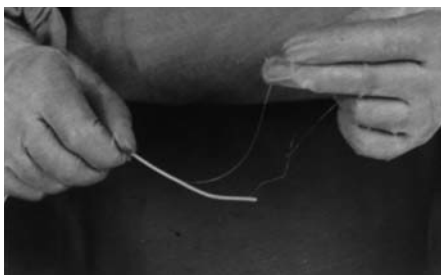In polypectomy snares, the loop end is already shaped. Insert the proximal loop end into the plastic probe from the distal end and push towards rear.

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Pull the end of the loop out of the plastic probe at the rear as far as possible, preferably arranging the probe in a circle. The free end of the loop must be long enough to extend below the locking screw.

Insert the end of the loop into the handle and lock it in place with the locking screw. Lock the probe connecting piece to the lower end of the handle.

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
Screw on thumb handle.
Check operation.

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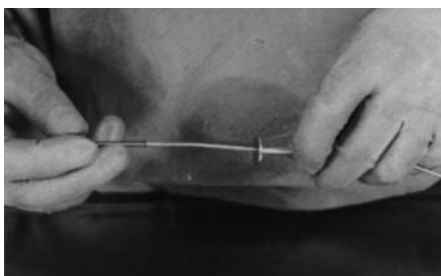
Flexible injection cannula

For preparation, pull the needle out of the flexible casing.

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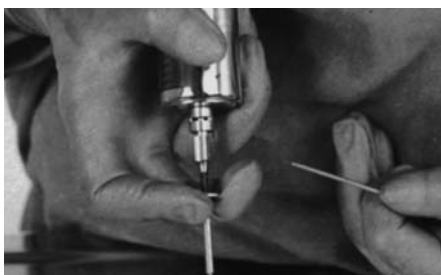Rinse the cladding, then immerse it in solution and clean in the ultrasound bath.

Rinse out the cannula using a syringe with disinfectant solution.

Afterwards, rinse with microbiologically pure/sterile water.

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Servicing and repair

Defective items of equipment must be serviced and repaired exclusively by persons authorized by us; all repair work must employ original parts.

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Limitation of liability

KARL STORZ GmbH & Co. KG shall be liable for failure or deterioration in the safe operation, operational reliability, and performance of this equipment only subject to the conditions that all assembly operations, system expansions, readjustments, modifications, or repairs to same have been performed by KARL STORZ, and that the instrument has been used in accordance with its operating instructions at all times.

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Warranty

For **two** years as of delivery to the end customer, we agree to replace the goods free of charge if proof can be provided of faulty materials or faulty workmanship. In doing so we cannot accept the cost of transportation or the risk of shipment. The warranty referred to in our Standard Conditions of Business shall apply.

Please fill out the attached warranty card on the last page and return it as soon as possible to:

KARL STORZ GmbH & Co. KG
Postfach 230
D-78503 Tuttlingen/Germany

Opening the equipment or performance of any repairs or modifications to the equipment by unauthorized persons shall release us from any liability for its performance. Any such opening, repair, or modification performed during the warranty period shall void all warranty.

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Directive compliance

This medical product bears the CE mark in accordance with the Medical Device Directive (MDD) 93/42/EEC. A code number after the CE mark indicates the responsible notified body.

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Repair and replacement program

KARL STORZ maintains a repair and replacement warehouse which is normally adequate to ensure prompt replacement of damaged telescopes and instruments. Under the repair and replacement plan, you receive an identical as-new instrument and are only charged the repair costs for the defective instrument. For telescopes, you receive a warranty of 1 year, and for instruments 6 months.

For **fiberscopes** and equipment, individual repair is necessary. To bridge the repair period, you will receive a unit on loan which you then return to KARL STORZ as soon as you receive the repaired unit.

In Germany you can refer repairs direct to

KARL STORZ GmbH & Co. KG
Abt. Reparaturservice
Dr. Karl-Storz-Strasse 34
D-78532 Tuttlingen/Germany

In other countries please contact your local KARL STORZ branch or authorized dealer.

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Important information

In order to prevent infection, it is strictly forbidden to ship contaminated medical devices.

All medical devices must be decontaminated on site, to avoid contact and aerogenous infections (among staff). We reserve the right to return contaminated instruments/equipment to the sender.

Repairs, changes or expansions which are not carried out by KARL STORZ or by experts authorized by KARL STORZ will invalidate all guarantee rights.

KARL STORZ gives no guarantee on the correct functioning of equipment or instruments which have been repaired by unauthorized third parties.

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Type

- 1 = Cleaner
2 = Disinfectant
3 = Suitable for ultrasound baths
A = Not suitable for flexible endoscopes (according to manufacturers' instructions)
(A) = Only for the flexible endoscopes of the new generation (year of manufacture > 2005)

Note: You will also find the current list of approved products on our Web site (www.karlstorz.com).

Note: When agents other than those named are used, KARL STORZ assumes no liability for any damage which may occur.

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*** Caution:** The use of CAUSTIC SODA can cause changes to the surface of aluminium parts (including coated parts), plastics and soldered connections, and may reduce the service life of the instrument set. For **flexible endoscopes**, only immerse the examination shaft, **not** the housing/handle part.

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Type

- 1 = Cleaner
2 = Disinfectant
3 = Suitable for ultrasound baths
4 = Neutralizer
5 = Rinse aid
A = Not suitable for flexible endoscopes
(according to manufacturers' instructions)

Note: You will also find the current list of approved products on our Web site (www.karlstorz.com).

Note: When agents other than those named are used, KARL STORZ assumes no liability for any damage which may occur.

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Machine preparation

The following instruments and telescopes produced by KARL STORZ must not be prepared in a cleaning and disinfecting machine: HAMOU® telescopes with ocular focusing, magnifiers with focusing, IMPERATOR drill handles, standard, tungsten carbide and diamond-tipped drills, prismatic light deflectors, and magnetic instrument rack. For motors and their handles the appropriate manual provides detailed information. Selection of a method for cleaning and disinfection must be discussed with the machine manufacturer and the manufacturer of the chemical reagents. Only special methods must be used which have been verified for this purpose. Regarding the microbiological effectiveness please ask the manufacturer. Thermal disinfection is preferable. The following agents for machine cleaning and disinfection have been approved regarding material compatibility:

Schülke & Mayr GmbH	Thermosept alka clean	1
	Thermosept DK	2
	Thermosept ED	2
	Thermosept ER	1
	Thermosept RKN-zym	1
Steris	Instru-Klenz	1, 3
Weigert, Dr. GmbH & Co.	neodisher DuoClean	1, A
	neodisher FA	1
	neodisher FA forte	1, 3
	neodisher Oxivario	1
	neodisher mediclean	1
	neodisher mediclean forte	1, 3
	neodisher mediklar	5
	neodisher medizym	1
	neodisher SeptoClean	1, 2, A
	neodisher Septo DA	2
	neodisher Septo DN	2

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Lubrication gels

Please note that the manufacturer's instructions for use must be closely observed when using the lubrication gels. We recommend the following lubrication gels:

Caution: Encrusted lubricant gel may damage the fiberscope. Remove any lubricant gel from the fiberscope surface by rinsing it immediately after use!

Note: Due to the special surface of intubation fiberscopes, KARL STORZ recommends the use of microbiologically pure/sterile water.

The following lubricant gels are safe for use with KARL STORZ fiberscopes and can be rinsed off.

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In the case of lubricants other than those listed above, KARL STORZ does not accept liability for any damage which may arise.

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WARRANTY

For a period of **two** years after delivery to the end-user, we shall replace free of charge equipment that can be proven to have defects in material or workmanship.

In doing so we cannot accept to bear the cost of transportation or the risk of shipment. The warranty referred to in our Standard Conditions of Business shall apply.

Complete warranty card upon purchase/delivery and return immediately to:

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Firma/Fecha:

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KARL STORZ GmbH & Co. KG
Postfach 230
78503 Tuttlingen/Alemania

To be filled out by instrument owner:

Return address/company stamp:

Please attach
sufficient
postage

Field of application:

Type of instrument:

Serial no.:

Purchase date:

Signature/Date:

REPLY CARD

KARL STORZ GmbH & Co. KG
Postfach 230
78503 Tuttlingen/Germany

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Absender/Firmenstempel:

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Geräte-Typ:

Serien-Nr.:

Kauf-Dat.:

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ANTWORTKARTE

KARL STORZ GmbH & Co. KG
Postfach 230
78503 Tuttlingen

KARL STORZ GmbH & Co. KG
Mittelstraße 8, 78532 Tuttlingen, Germany
Postfach 230, 78503 Tuttlingen, Germany
Phone: +49 (0)7461 708-0,
Fax: +49 (0)7461 708-105
E-Mail: info@karlstorz.de
Web: www.karlstorz.com

KARL STORZ Endoskope Berlin GmbH
Ohlauer Straße 43
10999 Berlin, Germany
Phone: +49 (0)30 30 69090
Fax: +49 (0)30 30 19452

KARL STORZ Endoscopy Canada Ltd.
2345 Argenta Road, Suite 100
Mississauga, ON, L5N 8K4, Canada
Phone: +1 905 816-8100, Fax: +1 905 858-0933
E-Mail: info@karlstorz.ca

KARL STORZ Endoscopy-America, Inc.
600 Corporate Pointe
Culver City, CA 90230-7600, USA
Phone: +1 310 338-8100, +1 800 421-0837
Fax: +1 310 410-5527
E-Mail: info@ksea.com

KARL STORZ Veterinary Endoscopy America, Inc.
175 Cremona Drive
Goleta, CA 93117, USA
Phone: +1 805 968-7776, Fax: +1 805 685-2588
E-Mail: info@karlstorvet.com

KARL STORZ Endoscopia Latino-America, Inc.
815 N. W. 57th Avenue, Suite 480
Miami, FL 33126-2042, USA
Phone: +1 305 262-8980, Fax: +1 305 262-8986
E-Mail: info@ksela.com

KARL STORZ Endoscopia
Miramar Trade Center
Edificio Jerusalem, Oficina 308 A,
La Habana, Cuba
Phone: +53 72041097, Fax: +53 72041098

KARL STORZ Endoscopia México S.A. de C.V.
Lago Constanza No 326
Col. Chapultepec Morales
D.F.C.P. 11520, México, Mexico
Phone: +52 5552 5056 07, Fax: +52 55 5545 0174

KARL STORZ Endoscopia Argentina S.A.
Cervino 4449 Piso 10°
1425 Buenos Aires C. F., Argentina
Phone: +54 11 4772 4545, Fax: +54 11 4772 4433
E-Mail: info@karlstorz.com.ar

KARL STORZ Endoskop Sverige AB
Storsåtragränd 14, 12739 Skärholmen, Sweden
Postal address: Po Box 8013, 14108 Kungens
Kurva, Sweden
Phone: +46 8 50 56 480, Fax: +46 8 50 56 4848
E-Mail: kundservice@karl-storz.se

KARL STORZ Endoskopi Danmark A/S
Skovlytoften 33
2840 Holte, Denmark
Phone: +45 45162600, Fax: +45 45162609
E-Mail: marketing@karlstorz.dk

KARL STORZ Endoscopy (UK) Ltd.
392 Edinburgh Avenue, Slough
Berkshire, SL1 4UF, Great Britain
Phone: +44 1753 503500, Fax: +44 1753 578124
E-Mail: customerservice@karlstorz-uk.com

KARL STORZ Endoscopia Nederland B. V.
Phone: +31 651 938 738, Fax: +31 135 302 231

KARL STORZ Endoscopy Belgium N. V.
Phone: +32 473 810 451, E-Mail: info@stopler.be

KARL STORZ Endoscopia France S. A.
12, rue Georges Guynemer
Quartier de l'Europe
78280 Guyancourt, France
Phone: +33 1 30484200, Fax: +33 1 30484201
E-Mail: marketing@karlstorz.fr

KARL STORZ Endoskop Austria GmbH
Landstraßer Hauptstr. 146/11/18
1030 Wien, Austria
Phone: +43 1 71 56 0470, Fax: +43 1 71 56 0479
E-Mail: storz-austria@karlstorz.at

KARL STORZ Endoscopia Italia S. r. l.
Via dell'Artigianato, 3
37135 Verona, Italy
Phone: +39 045 8222000, Fax: +39 045 8222001
E-Mail: info@karlstorz.it

KARL STORZ Endoscopia Ibérica S. A.
Parque Empresarial San Fernando
Edificio Francia – Planta Baja
28830 Madrid, Spain
Phone: +34 91 6771051, Fax: +34 91 6772981
E-Mail: marketing@karlstorz.es

KARL STORZ Endoskope Greece Ltd.*
Ipsilantou Str. 32
54248 Thessaloniki, Greece
Phone: +30 2310 304868 Fax: +30 2310 304862
*Repair & Service Subsidiary

KARL STORZ Industrial**
Gedik Is Merkezi B Blok
Kat 5, D 38-39, Bagdat Cad. No: 162
Maltepe Istanbul, Turkey
Phone: +90 216 442 9500, Fax: +90 216 442 9030
**Sales for Industrial Endoscopy

KARL STORZ Endoscopia ROMANIA srl
Str. Prof. Dr. Anton Colorian, nr. 74, Sector 4
041393 Bukarest, Romania
Phone: +40 (0)31 4250800, Fax: +40 (0)31 4250801
E-Mail: office@karlstorz.ro

TOV KARL STORZ Ukraine
18b Geroev Stalingrada avenu
04210 Kiev, Ukraine
Phone: +380 44 42668-14, -15, -19, -20
Fax: +380 44 42668-14

OOO KARL STORZ Endoscopy – WOSTOK
Derbenyevskaya nab. 7, building 4
115114 Moscow, Russia
Phone: +7 495 983 02 40, Fax: +7 495 983 02 41
E-Mail: kste-wostok@karlstorz.ru

KARL STORZ Endoskope
Regional Center for Endoscopy S.A.L.
St. Charles City Center, 5th Floor
Phoenicia Street, Mina Elhosh
2020 0908 Beirut, Lebanon
Phone: +961 1 368181, Fax: +961 1 365151

KARL STORZ Endoscopy South Africa, (Pty) Ltd.
8th Floor Convention Tower
Cnr Heerengracht & Coen Steytler, Foreshore
Cape Town 8001, South Africa
P.O.Box 6061, Roggebaai Cape Town,
8012, South Africa
Phone: +27 21 417 2600, Fax: +27 21 421 5103
E-Mail: info@karlstorz.co.za

TOO KARL STORZ ENDOSCOPY Kazakhstan
Khodjanova 17
050060 Almaty, Kazakhstan
Phone/Fax: +7 72 72 49 43 63, +7 72 72 49 41 00
e-mail: KarlStorzKazakhstan@itte.kz

KARL STORZ Endoscopy India Private Ltd.
C-126, Okhla Industrial Area, Phase-1
New Delhi 110020, India
Phone: +91 11 26815445-51,
Fax: +91 11 268129 86
E-Mail: karlstorz@vsnl.com

KARL STORZ GmbH & Co. KG
Resident Representative Office
80/33 (44/19) Dang Van Ngu,
F.10-Q. Phu Nhuan
Ho Chi Minh City, Vietnam
Phone: +848 991 8442, Fax: +848 844 0320

KARL STORZ Endoscopy (S) Pte. Ltd.
#05 – 08 San Centre, 171 Chin Swee Road
Singapore 0316, Singapore
Phone: +65 6376-1066, Fax: +65 6376-1068
E-Mail: karlstorz@pacific.net.sg

KARL STORZ Endoscopy Australia Pty. Ltd.
174 Parramatta Rd.
Camperdown, N.S.W. 2050, Australia
Phone: +61 02 8594 9150, Fax: +61 02 8594 9199
E-Mail: info@karlstorz.au

KARL STORZ Endoscopy China Ltd.
Hong Kong Representative Office
Unit 1601, Chinachem Exchange Square
1 Hoi Wan Street, Quarry Bay,
Hong Kong, People's Republic of China
Phone: +8 52 28 65 2411, Fax: +8 52 28 65 4114
E-Mail: inquiry@karlstorz.com.hk

KARL STORZ Endoscopy China Ltd.
Beijing Representative Office
Room 610, China Life Tower
No. 6, Chaowai Street
Beijing, 100020, People's Republic of China
Phone: +86 10 8525 3725, Fax: +86 10 8525 3728
E-Mail: inquiry@karlstorz.com.hk

KARL STORZ Endoscopy China Ltd.
Unit 3901-3904, Tower 1 Grand Gateway
No. 1 Hong Qiao Road
Shanghai, 200030, People's Republic of China
Phone: +86 21 6113-1188, Fax: +86 21 6113-1199
E-Mail: inquiry@karlstorz.com.hk

KARL STORZ Endoscopy China Ltd.
Chengdu Representative Office
F-5, 24/F., Chuanxing Mansion,
No. 18 Renming Road South
Chengdu, Sichuan, 610016, People's Rep. of China
Phone: +86 28 8620-0175, Fax: +86 28 8620-0177
E-Mail: inquiry@karlstorz.com.hk

KARL STORZ Endoscopy China Ltd.
Guangzhou Representative Office
Room 1119-20, Dongshan Plaza
69 Xianlie Road Middle, Dongshan District,
Guangzhou, Guangdong, 510095,
People's Rep. of China
Phone: +86 20 8732-1281, Fax: +86 20 8732-1286
E-Mail: inquiry@karlstorz.com.hk

KARL STORZ Endoscopy Asia Marketing Pte Ltd
3791 Jalan Bukit Merah
06-11 e-Centre @ Redhill
Singapore 159471, Singapore
Tel. No. +65 63761066, Fax. No. +65 63761068
Email : infoasia@karlstorz.com.sg
Email : serviceasia@karlstorz.com.sg

KARL STORZ Endoscopy Singapore Sales Pte Ltd
3791 Jalan Bukit Merah
10-18 e-Centre @ Redhill
Singapore 159471
Tel. No. +65 65325548, Fax No. +65 65323832
Email : infosg@karlstorz.com.sg
Email : servicesg@karlstorz.com.sg

KARL STORZ Endoscopy Japan K. K.
Bois Hongo Building 6Fl, 3-42-5 Hongo
Bunkyo-ku, Tokyo 113-0033, Japan
Phone: +81 3 5802-3966, Fax: +81 3 5802-3988
E-Mail: info@karlstorz.jp

www.karlstorz.com

KARL STORZ GmbH & Co. KG

Mittelstraße 8
78532 Tuttlingen

Postfach 230
78503 Tuttlingen
Germany

Telefon: +49 (0)7461 708-0
Telefax: +49 (0)7461 708-105
E-Mail: info@karlstorz.de
Web: www.karlstorz.com